

## Recombinant Cluster Of Differentiation 14 (CD14)

Catalog No.: TP01478

50µg

### Sequence Information

**Species:** Human

**Gene ID:**929

**Swiss Prot:**P08571

**Synonyms:**Monocyte differentiation antigen

CD14; Myeloid cell-specific leucine-rich glycoprotein

**Residues:**Asp125~Leu288

DLKITGTMPPLEATGLALSSLRLRNVSATGRSWLAELQQWLKPKLVLSIA

QAHSFAFSCEQVRAFPALTSLDLSDNPGLGERGLMAALCPHKFPAIQNLALRNT

GMETPTGVCAALAAAGVQPHSLDLSHNSLRATVNPSPAPRCMWSSALNSLNLSFA

GL

### Product Information

**Source:** Prokaryotic expression.

**Host:** *E. coli*

**Tags:**N-terminal His Tag

**Subcellular Location:**Secreted.

**Purity:** >97%

**Traits:** Freeze-dried powder

**Buffer formulation:**PBS, pH7.4, containing 0.1% SKL, 5% Trehalose.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 8.5

**Predicted Molecular Mass:** 18.6kDa

**Accurate Molecular Mass:** 18kDa as determined by SDS-PAGE reducing conditions.

### [ USAGE ]

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

### [ STORAGE AND STABILITY ]


**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.


Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

**[ IDENTIFICATION ]**

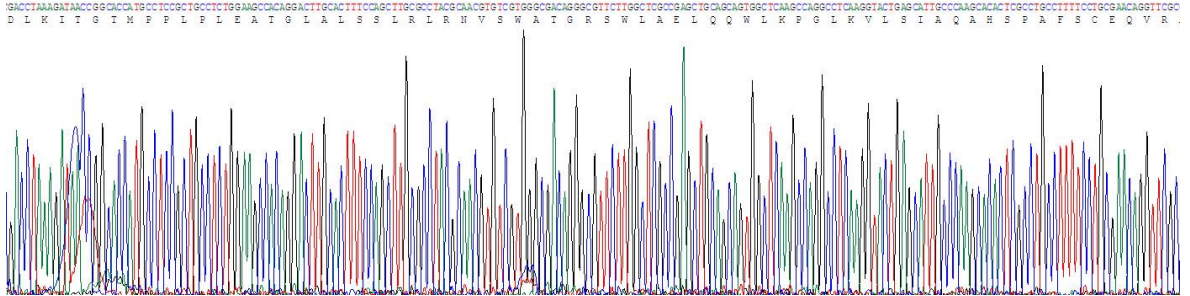

**Figure 1. Gene Sequencing (Extract)**

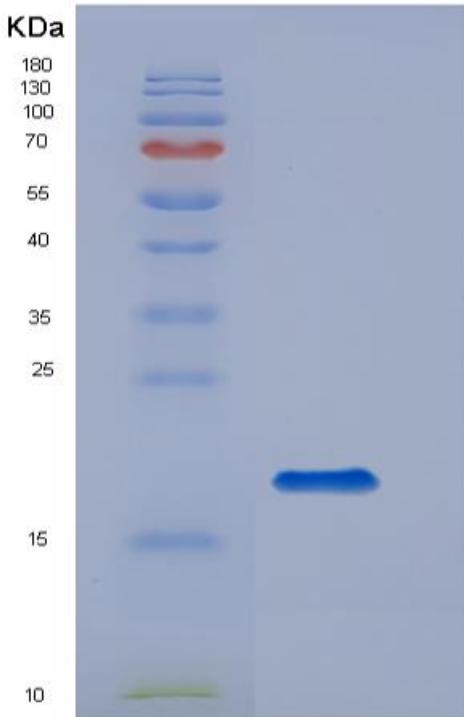

**Figure 2. SDS-PAGE**

**[ IMPORTANT NOTE ]**



The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.