

## Polyclonal Antibody to Bcl2 Associated X Protein (Bax)

PA03835      100ul

### **[ PROPERTIES ]**

**Source:** Polyclonal antibody preparation

**Host:** Rabbit

**Purification:** Recombinant Protein affinity chromatography.

**Traits:** Liquid

Concentration: 500ug/ml

Ig Isotype:IgG

**UOM:** 100μl

**Cross Reactivity:** Human, Mouse, Rat

**Applications:** WB, IHC, IP;Elisa.

### **[ IMMUNOGEN ]**

Immunogen: Recombinant protein of human Bax.

### **[ APPLICATIONS ]**

**Western blotting:** 0.5-2μg/mL;

**Immunohistochemistry:** 5-20μg/mL;

**Immunocytochemistry:** 5-20μg/mL;

Optimal working dilutions must be determined by end user.

### **[ FORMULATION ]**

**Form & Buffer:** Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05%

Proclin-300, 50% glycerol.

### **[ STORAGE AND STABILITY ]**


**Storage:** Avoid repeated freeze/thaw cycles.

Store at 4°C for frequent use.

Aliquot and store at -20°C for 12 months.


**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

### [ IDENTIFICATION ]

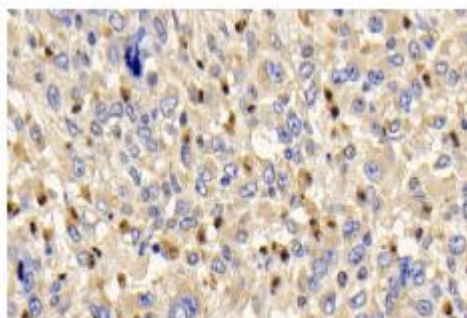

DAB staining on IHC-P;

Samples: Human Liver cancer Tissue;

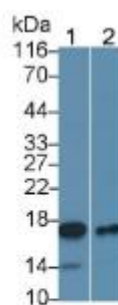

Western Blot; Sample: Lane1: Human Lung  
lysate; Lane2: Porcine Lung lysate Primary

### [ IMPORTANT NOTE ]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.

---