

Recombinant Interleukin 2 (IL2)

Catalog No.: TP04711

100µg

Sequence Information

Species: Human**Gene ID:**3558**Swiss Prot:**P60568**Synonyms:**TCGF**Residues:**Ala21~Thr153

```
          APTSSSTKKT QLQLEHLLLD LQMILNGINN
YKNPKLTRML TFKFYMPKKA TELKHLQCLE EELKPLEEVL NLAQSKNFHL
RPRDLISNIN VIVLELKGSE TTFMCEYADE TATIVEFLNR WITFCQSIIS
TLT
```

Product Information

Source: Prokaryotic expression..**Host:** E. coli**Tags:** N-terminal His-Tag**Subcellular Location:** Secreted.**Purity:** >90%**Traits:** Freeze-dried powder**Buffer formulation:** 20mM Tris, 150mM NaCl, pH8.0, 5% Trehalose.**Original Concentration:** 200µg/mL**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 7.1**Predicted Molecular Mass:** 15.0kDa**Accurate Molecular Mass:** 15kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one week.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[**IDENTIFICATION**]

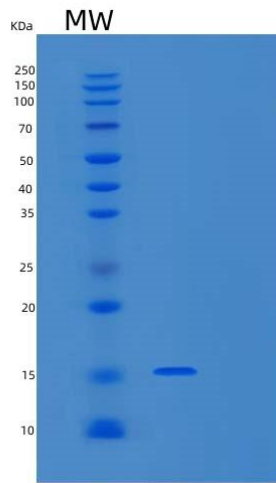


Figure 2. SDS-PAGE

[**IMPORTANT NOTE**]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.