

Eukaryotic A Disintegrin And Metalloprotease 10 (ADAM10)

Catalog No.: TP01058

50µg

Sequence Information

Species: Human

Gene ID:102

Swiss Prot:O14672

Synonyms:Disintegrin and metalloproteinase domain-containing protein 10 isoform 1 , ADAM10, AD10, AD18, CD156c, HsT18717, kuz, MADM, RAK

Residues:Thr214~Glu672

TTSAEKNTCQLYIQTDHLFFKYVGTREAVIAQISSHVKAIDTIYQTTDFSGIRN
ISFMVKRIRINTTADEKDPTNPF RFPNIGVEKFLELNSEQNHDDYCLAYVFTDR
DFDDGVLGLAWVGAPSGSSGGICEKSKLYSDGKKKSLNTGIITVQNYGSHVPPK
VSHITFAHEVGHNF GSPHDSGTECTPGESKNLGQKENGNYIMYARATSGDKLNN
NKFSLCSIRNISQVLEKRRNNCFVESGQPICGNGMVEQGEEDCGYSDQCKDEC
CFDANQPEGRKCKLKPGKQCSPSQGPCCTAQCAFKSKSEKCRDDSDCAREGICN
GFTALCPASDPKPNFTDCNRHTQVCINGQCAGSICEKYGLEECTCASSDGKDDK
ELCHVCCMKKMDPSTCASTGSVQWSRHFSGRITITLQPGSPCNDFRGYCDVFMRC
RLVDADGPLARLKKAI FSPELYENIAE

Product Information

Source: Eukaryotic expression.

Host: CHO Cell

Tags: N-terminal His-Tag

Subcellular Location: Membrane.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.2

Predicted Molecular Mass: 52.2kDa

Accurate Molecular Mass: 52kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

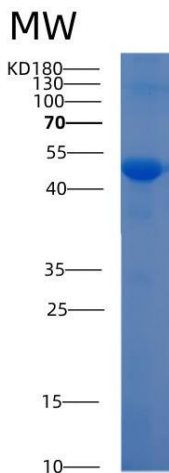


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.