



Recombinant Human V-Ki-Ras2 Kirsten Rat Sarcoma Viral Oncogene Homolog (KRAS)

Catalog No.: TP02762 50µg

Sequence Information

Species: Human

Gene ID:3845

Swiss Prot:P01116

Synonyms:GTPase KRas isoform a precursor,
C-K-RAS, K-RAS2A, K-RAS2B,
K-RAS4A, K-RAS4B, KI-RAS, KRAS1,
KRAS2, NS, NS3, RASK2

Residues:Met1~Met189

MTEYKLVVVGAGGVGKSALTIQLIQNHFVDEYDPTIEDSYRKQVVIDGETCLLD
ILDTAGQEEYSAMRDQYMRTGEGFLCVFAINNTKSFEDIHHYREQIKRVKDS
E
VPMVLVGNKCDLPSRTVDTKQAQDLARSYGIPFIETSAKTRQRVEDAFYTLVRE
IRQYRLKKISKEEKTGCVKIKKCIIM

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Membrane

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.3

Predicted Molecular Mass: 25.4kDa

Accurate Molecular Mass: 25kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

TCCTGGGACCCATAAGCCTAGCATGCTGGTGGACGAAATGGGTGGGGATCCATGACTGAAATAACTTTGTGGTGGTGGACCTGGTGGGCTAGCCAGAGTGGCTTTGACGATACAGCTAATTCGAGATCATTTTGTGGGCAATATGATCGAACATAGAGGATTCCTACAGGAAAGCAAGTAGTAATTTGATGGAGAACTGT

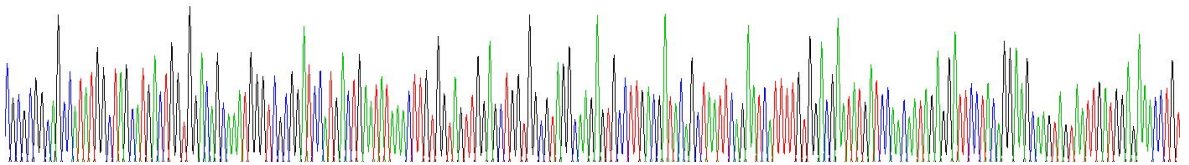


Figure 1. Gene Sequencing (Extract)

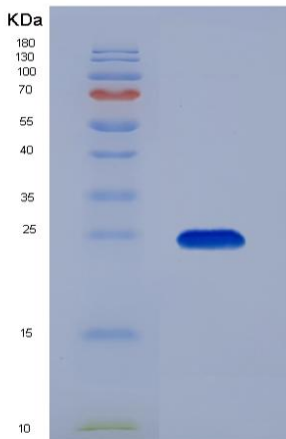


Figure 2. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.