

Recombinant Human Cluster Of Differentiation 4 (CD4)

Catalog No.: TP01513

50µg

Sequence Information

Species: Human

Gene ID:920

Swiss Prot:P01730

Synonyms:T-cell surface glycoprotein CD4
isoform 1, CD4mut

Residues:Lys26~Trp390

KKVVLGKKGDTVELTCTASQKKSIFHWKNSNQIKILGNQGSFLTKGPSKLNDR
ADSRRLWDQGNFPLIIKNLKIEDSDTYICEVEDQKEEVQLLVFGLTANS DTHL
LQGQSLTLTLESPPGSSPSVQCRSPRGKNIQGGKTL SVSQLELQDSGTWTCTVL
QNQKKVEFKIDIVVLAFAQASSIVYKKEGEQVEFSFPLAFTVEKLTGSGELWWQ
AERASSKSWITFDLKNKEVSVKRVTQDPK LQMGKKLPLHLTLPQALPQYAGSG
NLTLALEAKTGK LHQEVNLVVMRATQLQKNLTCEVWGPTSPKLM LSLKLENKEA
KVS KREKAVWVLNPEAGMWQCLLSDSGQV LLESNIKVLPTW

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Membrane

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 9.2

Predicted Molecular Mass: 41.7kDa

Accurate Molecular Mass: 44kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

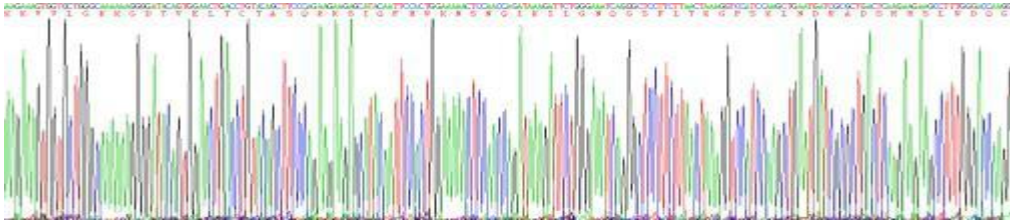


Figure 1. Gene Sequencing (Extract)

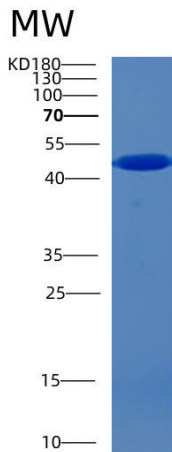


Figure 2. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.