

Eukaryotic Alpha-1-Microglobulin (a1M)

Catalog No.: TP06475

50µg

Sequence Information

Species: Human

Gene ID:259

Swiss Prot:P02760

Synonyms:A1M;EDC1;HCP;HI30;IATIL;ITI;ITIL;ITILC;UTI

Residues:Gly20~Asn352

GPVPTPPDNIQVQENFNISRIYGKWYNLAIGSTCPWLKKIMDRMTVSTLVLGEG
ATEAEISMTSTRWRKGVCEETSGAYEKTDGKFLYHKSQWITMESYVVHTNY
DEYAIFLTKKFSRHHGPTITAKLYGRAPQLRETLQDFRVVAQGVGIPEDSIFT
MADRGEVCPGEQEPEPILIPRVRAVLPQEEEGSGGGQLVTEVTKKEDSCQLGY
SAGPCMGMTSRFYNGTSMACETFQYGGCMGNGNMFVTEKECLQTCRTVAACNL
PIVRGPCRAFIQLWAFDAVKGKCVLFPYGGCQGNKGFYSEKECREYCGVPGDG
DEELLRFSN

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: N-terminal His Tag and C-terminal Fc Region of Human IgG1

Subcellular Location: Secreted.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 1000µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.8

Predicted Molecular Mass: 64.3kDa

Accurate Molecular Mass: 54kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

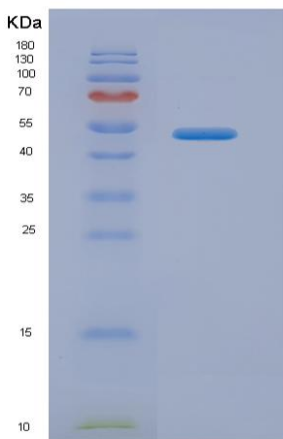


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.