



## Recombinant Mouse Annexin A1 (ANXA1)

Catalog No.: TP01178 1mg

### Sequence Information

**Species:** Mouse

**Gene ID:**16952

**Swiss Prot:**P10107

**Synonyms:**Annexin A1, Anx-1, Anx-A1,  
C430014K04Rik, Lpc-1, Lpc1

**Residues:**Glu6~Asn346

EFLKQARFLENQEYVQAVKSYKGGPGSAVSPYPSFNVSSDVAALHKAIMVKGVD EATIIDILTKRTNAQRQQI  
KAAYLQENGKPLDEVLRKALTGHLEEVVLAMLKTPAQFDADELRGAMKGLGTDEDTLIEILTTRSNEQIREINRV  
YREELKRD LAKDITSDTSGDFRKALLALAKGDRCQDLSVNQDLADTDARALYEAGERRKGT DVNVFTTILTSR SF  
PHLRRVFQNYGKYSQHDMNKALDLELKGDI EKCLTTIVKCATSTPAFFAEKLYEAMKGAGTRHKALIRIMVSRSE  
IDMNEIKVFYQKKYGISLCQA ILDETKGDYEKILVALCGGN

### Product Information

**Source:** Prokaryotic expression.

**Host:** *E. coli*

**Tags:** N-terminal His-Tag

**Subcellular Location:** Membrane, Nucleus, Secreted, Cytoplasm.

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS pH7.2, 0.01% sarcosyl, 5%Trehalose .

**Original Concentration:** 1mg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 6.2

**Predicted Molecular Mass:** 44.9kDa

**Accurate Molecular Mass:** 44kDa as determined by SDS-PAGE reducing conditions.

### [ USAGE ]

Reconstitute in PBS pH7.2 or ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

### [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

GGGATCGAATTCCTCAAGGCGGCGGTTTCTTGAAATCAAGACGGGATATGTTCAAGCTGAATCATCAAGGAGGCTCTGGGTCAGGAGTGAAGGCGGCTTCCCTTCAATGATATCTGGATGCTCTGGCTTGGACAAAGCTATCATGGTTAAGGCTGATGAGCAACATCATTCATCTTAAAGGAGGAG

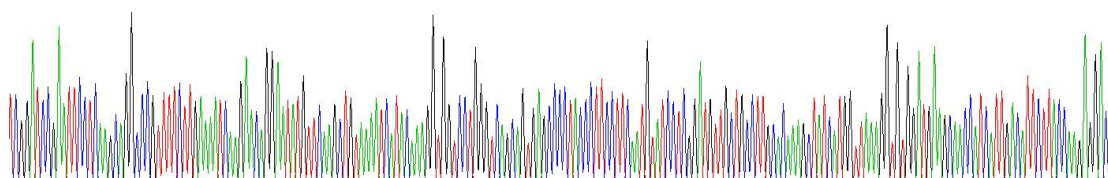


Figure 1. Gene Sequencing (Extract)

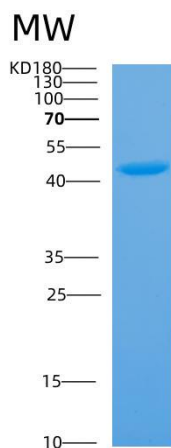


Figure 2. SDS-PAGE

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.