

Recombinant Human Non Metastatic Cells 1, Protein NM23A Expressed In (NME1)

Catalog No.: TP05739 100µg

Sequence Information

Species: Human

Gene ID:4830

Swiss Prot:P15531

Synonyms:Nucleoside Diphosphate Kinase A;
NDK A; NDP Kinase A; Granzyme
A-Activated DNase; GAAD; Metastasis
Inhibition Factor nm23; Tumor
Metastatic Process-Associated Protein;
nm23-H1; NME1; NDPKA; NM23

Residues:Met 1-Glu152

MANCERTFIAIKPDGVQRGLVGEIIKRFEQKGFRLVGLKFMQASEDLLKEHYVD

LKDRPFFAGLVKYMHS GPVAMVWEGLNVVKTGRVMLGETNPADSKPGTIRGDF

CIQVGRNIIHGSDSVESAEKEIGLWFHPEELVDYTSCAQNWIYE

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Cytoplasm. Nucleus

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.6

Predicted Molecular Mass: 20.8kDa

Accurate Molecular Mass: 22kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

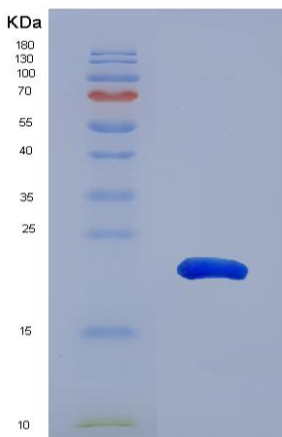


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.