

Recombinant Synaptosomal Associated Protein 25kDa (SNAP25)

Catalog No.: TP03977

100µg

Sequence Information

Species: Human

Gene ID:6616

Swiss Prot:P60880

Synonyms:SNAP25, SNAP, SNAP25,
Synaptosomal-associated protein, His-
tag, SNAP-25, Synaptosomal-associated
25 kDa protein, Super protein, SUP,
Synaptosomal-associated protein 25,
Synaptosomal-associated protein 25
isoform SNAP25B bA416N4.2

Residues:Met1~Gly206

MAEDADMRNELEEMQRRADQLADESLESTRRMLQLVEESKDAGIRTLVMLDEQG
EQLERIEEGMDQINKDMKEAEKNLTDLGKFCGLCVPCNKLKSSDAYKKAWGNN
QDGVVASQPARVVDEREQMAISGGFIRRVTDARENEMDENLEQVSGIIGNLRH
MALDMGNEIDTQNRQIDRIMEKADSNKTRIDEANQRATKMLGSG

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His Tag

Subcellular Location:Nucleus.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.2

Predicted Molecular Mass: 19.8kDa

Accurate Molecular Mass: 20kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

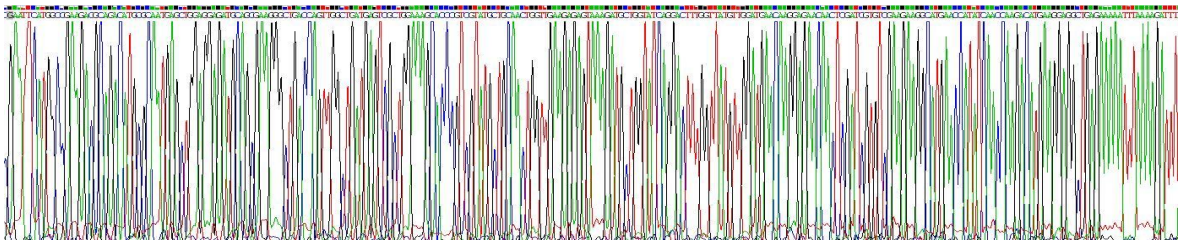


Figure 1. Gene Sequencing (Extract)

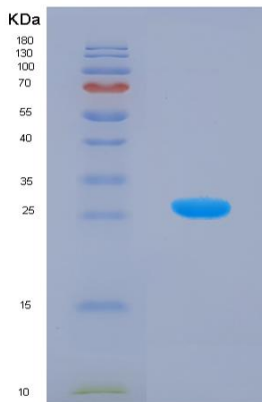


Figure 2. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.