

## Recombinant Activity Regulated Cytoskeleton Associated Protein (ARC)

Catalog No.: TP01204 100µg

### Sequence Information

**Species:** Human

**Gene ID:**23237

**Swiss Prot:**Q7LC44

**Synonyms:**Activity-regulated

cytoskeleton-associated protein, Arg3.1,

KIAA0278

**Residues:**Cys96~Gln356

CRCQETIANLERWVKREMHVWREVFYRLERWADRLESTGGKYPVGSESARHTVS  
VGVGGPESYCHEADGYDVTSPYAITPPPAAGELPGQEPAEAQQYQPWVPGEDG  
QPSPGVDTQIFEDPREFLSHLEEYLRQVGGSEEWLSQIQNHMNGPAKKWWEFK  
QGSVKNWVEFKKEFLQYSEGTL SREAIQRELDLPQKQGEPLDQFLWRKRDLYQT  
LYVDADEEEIIQYVVGTLPKLRFLRHPLPKTLEQLIQRGMEVQ

### Product Information

**Source:** Recombinant expression.

**Host:** *E.coli*

**Tags:** N-terminal His Tag

**Subcellular Location:** Chromosome.

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 6.1

**Predicted Molecular Mass:** 36.6kDa

**Accurate Molecular Mass:** 40kDa as determined by SDS-PAGE reducing conditions.

### [ USAGE ]

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

### [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

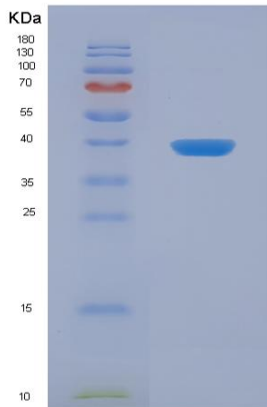


Figure 1. SDS-PAGE

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.