

Recombinant Annexin A10 (ANXA10)

Catalog No.: TP01179

100µg

Sequence Information

Species: Human

Gene ID:11199

Swiss Prot:Q9UJ72

Synonyms:Annexin A10, ANX14.

Residues:Met1~Tyr324

MFCGDYVQGTIFPAPNFNPIMDAQMLGGALQGFDCDKDMLINILTQRCNAQRMM
IAEAYQSMYGRDLIGDMREQLSDHFKDVMAGLMYPPPLYDAHELWHAMKGVGTD
ENCLIEILASRTNGEIFQMREAYCLQYSNNLQEDIYSETSGHFRDTLMNLVQGT
REEGYTDPAMAAQDAMVLWEACQQKTGEHKTMLQMILCNKSYQQLRLVVFQEFQN
ISGQDMVDAINECYDGYFQELLVAIVLVCVRDKPAYFAYRLYSIAIHDFGFHNKTV
IRIL IARSEIDLLTIRKRYKERYGKSLFHDIRNFASGHYKALLAICAGDAEDY

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His Tag

Subcellular Location: Mitochondrion

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.13

Predicted Molecular Mass: 41.0kDa

Accurate Molecular Mass: 38kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

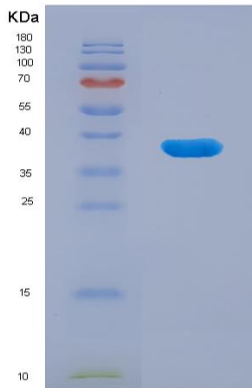


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.