

## Recombinant Human Zinc- $\alpha$ -2-Glycoprotein/AZGP1/ZAG Protein(C-6His)

Catalog No.: TP04882

50 $\mu$ g

### Sequence Information

**Species:** Human

**Gene ID:**563

**Swiss Prot:**P25311

**Synonyms:**Zinc-Alpha-2-Glycoprotein;

Zn-Alpha-2-GP;

Zn-Alpha-2-Glycoprotein; AZGP1; ZAG;

ZNGP1

**Residues:**Gln21-Ser298

QENQDGRYSLTYIYTGLSKHVEDVPAFQALGSLNDLQFFRYNSKDRKSQPMGLW  
RQVEGMEDWKQDSQLQKAREDFMETLKDIVEYYNDSNGSHVLQGRFGCEIENN  
RSSGAFWKYYYDGDYIEFNKEIPAWVPFDPAAQITKQKWEAEPVVVQRAKAYL  
EEECPATLRKYLKYSKNILDRQDPPSVVVTSHQAPGEKKLKLCLAYDFYPGKID  
VHWTRAGEVQPELRGBVLHNGNGTYQSWVVAVPPQDTAPYSCHVQHSSLAQP  
LVVPWEAS

### Product Information

**Source:** Eukaryotic expression.

**Host:** 293F cell

**Tags:** C-terminal His Tag

**Subcellular Location:** Secreted .

**Purity:** >95%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

**Original Concentration:** 1000 $\mu$ g/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 5.7

**Predicted Molecular Mass:** 33.3kDa

**Accurate Molecular Mass:** 42kDa as determined by SDS-PAGE reducing conditions.

### [ USAGE ]

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

## [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

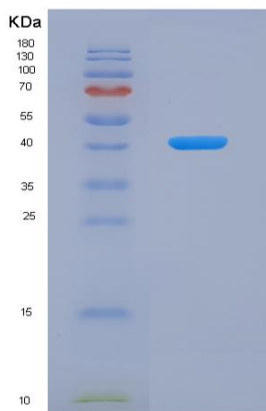


Figure 1. SDS-PAGE

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.