

Recombinant Human Visceral Adipose Tissue Derived Serine Protease Inhibitor (Vaspin)

Catalog No.: TP07034 50µg

Sequence Information

Species: Human

Gene ID:145264

Swiss Prot:Q8IW75

Synonyms:SERPINA12; OL-64; Visceral Adipose-Specific SERPIN; Serine Proteinase Inhibitor,Clade A Member 12; Visceral adipose-specific serpin

Residues:Leu20-Lys414

LLKPSFSRPNYKALSEVQGQWQRMAAKELARQNMDLGFKLLKKLAFYNPGRNIF
LSPLSISTAFSMLCLGAQDSTLDEIKQGFNFRKMPEKDLHEGFHYIHELTKQT
QDLKLSIGNTLFIDQRLQPQRKFLEDAKNFYSAETILTNFQNLMAQKQINDFI
SQKTHGKINLLIENIDPGTVMLLANYIFFRARWKHEFDPNVTKEEDFFLEKNSS
VKVPMFMRSGIYQVGYDDKLSCTILEIPYQKNITAIIFILPDEGKCLKHLEKGLQV
DTFSRWKTLSSRRVVDVSVPRLHMTGTFDLKKTLSYIGVSKIFEEHGDLTKIAP
HRSLKVGEAVHKAELKMDERGTEGAAGTGAQTLPMETPLVVKIDKPYLLLIYSE
KIPSVLFLGKIVNPIGK

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: C-terminal His-Tag

Subcellular Location: Cell membrane; Lipid-anchor, GPI-anch.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 1000µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 9.8

Predicted Molecular Mass: 54kDa

Accurate Molecular Mass: 54kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

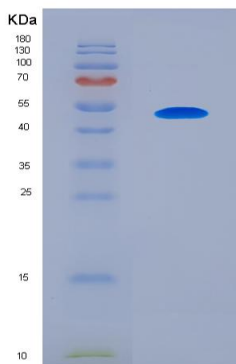


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.