



Recombinant Ferroptosis suppressor protein 1 (FSP1)

Catalog No.: TP11033

50µg

Sequence Information

Species: Rat

Gene ID:361843

Swiss Prot:A0A8I5ZQB7

Synonyms:AIFM2; p53-responsive gene 3

Residues:Met1-Pro403

MGTPRPRASVPGKRADRYRIFSQPSSAKFKMGSQVSVDTGAVHVVIVGGGFGGIAAASQLQALNVPFMLVDMKDSFHHNVAALRASVESGFAKKTFFIS
YSVTFKDNFRQGVIGIDLKNRMVLLLEGGEALPFSHLILATGSTGPFPGKFNEVSCQQAIIQAYEDMVKQIQRSQFIVVVGGSAGVEMAAEIKTEYP
EKEVTLIHSRVPLADKELLPCVRQEVKEILLRKGVQLLLSERVSNLEELPRNEYREYIKVETDKGTEVATNMVIVCNGIKINSSAYRSAFESRLASNG
ALKVNEFLQVEGYSNIYAIGDCADIKPEKMAHYAGLHANIAVANIVNSMKQRPLKAYKPGALTFLLSMGRNDGVGQISGFYVGRMLMVRLLAKSRDLLIS
TSWKTMRQSPP

Product Information

Source: Prokaryotic expression.

Host: *E. coli*

Tags: N-terminal His

Subcellular Location: Membrane.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% skl,5%Trehalose.

Original Concentration: 100µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 9.0

Predicted Molecular Mass: 45.2kkDa

Accurate Molecular Mass: 46kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

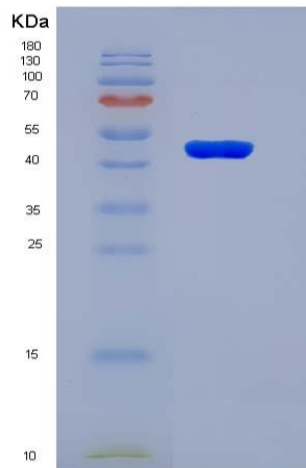


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.