

Recombinant Human Aromatase (CYP19A1)

Catalog No.: TP10261

100µg

Sequence Information

Species: Human

Gene ID:1588

Swiss Prot:P11511

Synonyms:CYPXIX;Cytochrome

P-450AROM;Cytochrome P450 19A11

Publication;Estrogen synthase

Residues:Met1-His503

MVLEMLNPIHYNITSIVPEAMPAATMPVLLLTGLFLLVWNYEGTSSIPGPGYCM
GIGPLISHGRFLWMGIGSACNYNRYVGEFMRVWISGEETLIISKSSSMFHIMK
HNHYSSRFGSKLGLQCIGMHEKGIIFNNNPELWKTRPFFMKALSGPGLVRMVT
VCAESLKTHLDRLEEVTNESGYVDVLTLLRRVMLDTSNTLFLRIPLDESAIVK
IQGYFDAWQALLIKPDIFFKISWLYKKYEKSVKDLKDAIEVLIAEKRRRISTEE
KLEECMDFATELILAEKRGDLTRENVNQCILEMLIAAPDTMSVSLFFMLFLIAK
HPNVEEAIKEIQTVIGERDIKIDDIQKLKVMENFIYESMRYQPVVDLVMRKAL
EDDVIDGYPVKKGTNIILNIGRMHRLEFFPKPNEFTLENFAKNVPYRYFQPFQGF
GPRGCAGKYIAMVMMKAILVTLLRRFHVKTLQGQCVEISQKIHDLSLHPDETKN
MLEMIFTPRNSDRCLEH

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Endoplasmic reticulum membrane.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 7.6

Predicted Molecular Mass: 62.0kDa

Accurate Molecular Mass: 62kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

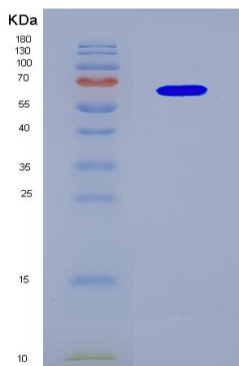


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.