



Recombinant Ubiquitin carboxyl-terminal hydrolase 5 (USP5)

Catalog No.: TP09938

100µg

Sequence Information

Species: Human

Gene ID:8078

Swiss Prot:P45974

Synonyms: ISOT

Residues:Met1-Ser858

MAELSEEALLSVLPTIRVPKAGDRVHKDECAFSFDTPSEGGLYICMNTFLGFGKQYVERHFNKTGQRVYLHLRRTRRPKEEDPATGTGDPPrKKPTRLAIG
VEGGFDLSEEFELDEDVKIVILPDYLEIARDGLGGLPDIVRDRVTSAVEALLSADSASRKQEVQAWDGEVRQVSKHAFSLKQLDNPARIPPCGWKCSKCDM
RENWLNLTDGSILCGRRYFDGSGGNNHAVEHYRETGYPLAVKLGITITPDGADVYSYDEDDMVLDPSLAEHLSHFGIDMLKMQKTDKTMTELEIDMNQRIGE
WELIQESGVPLKPLFGPGYTGIRNLGNSCYLNSVVQVLFSDIPDFQRKYVDKLEKIFQNAPTDPTQDFSTQVAKLGHGLLSGEYSKVPVPSGDGERVPEQKEV
QDGIAPRMFKALIGKGHPEFSTNRQQDAQEFFLHLINMVERNCRSSENPNVEFRFLVEEKIKCLATEKVYKTQRVDYIMQLPVPMDAALNKEELLEYYEKKR
QAEEEKMALPELVRAQVPFSSCLEAYGAPEQVDDFWSTALQAKSVAVKTTTFASFPDYLVIQIKKFTGLDWPKKLDVSIEMPEELDISQLRGTGLQPGEE
ELPDIAPPLVTPDEPKGSLGFYGNEDSFCSPHFSSPTSPMLDESVIIQLVEMGFPMDACRKAVYYTGNSGAEAMNWMVSHMDDPDFANPLILPGSSGPG
STAAADPPPEDCVTTIVSMGFSRDQALKALRATNNSLERAVDWIFSHIDDLDAEAAMDISEGRSAADSISESPVGPVKVRDGPQYQLFAFISHMGTSTMC
GHYVCHIKKEGRWVIYNDQKVCASEKPPKDLGYIYFYQRVAS

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-his and SUMO tag

Subcellular Location: Secreted.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 4.92

Predicted Molecular Mass: 96.0kDa

Accurate Molecular Mass: 96kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

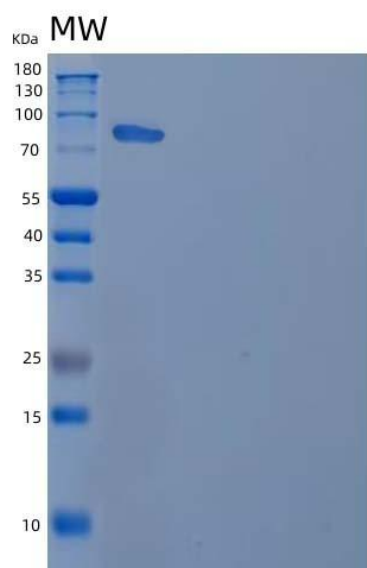


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.