

Recombinant Thyroglobulin (TG)

Catalog No.: TP11283

100µg

Sequence Information

Species: Pig

Gene ID: 100156471

Swiss Prot:F1RRV3

Synonyms:AITD3; TGN

Residues:Pro728-Cys923

PTPCQLQAEQAFGLGTVRGLISNPSEPPVLSSIIYPQCSASGQWRRVQCDGPPEQ

AFEWYERWGAQSRSGQELTPAELLMKIMSYREAASGSFRLFIQNLYEAGQQGIF

PGLARYSSLQDVPLAVLEGNTQATGNILLEPYLFWQILNGQLPRYPGPYSDFS

APLAHLDLRSCWCVDEAGRKLEGTQTEPSKVPAC

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Secreted.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 4.6

Predicted Molecular Mass: 25.4kDa

Accurate Molecular Mass: 26kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was

determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

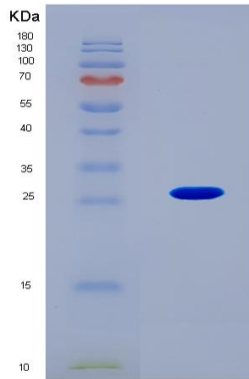


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.