

Recombinant Flavin-containing monooxygenase 3 (FMO3)

Catalog No.: TP10366

100µg

Sequence Information

Species: Human

Gene ID:2328

Swiss Prot:P31513

Synonyms:

Residues:Met1-Thr532

MGKKVAIIGAGVSGLASIRSCLEEGLEPTCFEKSNDIGGLWKFSDHAE EGRASI
YKSVFSNSSKEMMCFDPDFPDDFPNFMHNSKIQEYIIAFAKEKNLLKYIQFKT
FVSSVNKHPDFATTGQWDVTTTERDGKKESAVFDAMVCSGHHVYPNLPKESFPG
LNHFKGKCFHSRDYKEPGVFNGKRVLVVGLGNSGCDIATELSRTAEQVMISSRS
GSWMSRVWDNGYPWDMLLVTRFGTFLKNNLPTAISDWLYVKQMNARFKHENYG
LMPLNGVLRKEPVFNDELPASILCGIVSVKPNVKEFTETSAIFEDGTIFEGIDC
VIFATGYSFAYPFLDESIIKSRNNEIILFKGVFPPLLEKSTIAVIGFVQSLGAA
IPTVDLQSRWAAQVIKGTCTLPSMEDMMNDINEKMEKKRKWFGKSETIQTDYIV
YMDELSSEFIGAKPNIPWLFLTDPKLAMEVYFGPCSPYQFRLVGPQWPGARNAI
LTQWDRSLKPMQTRVVGRLLQKPCFFHWWLKLFAIPILLIAVFLVLT

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His and SUMO Tag

Subcellular Location: Microsome membrane.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 7.9

Predicted Molecular Mass: 76.0kDa

Accurate Molecular Mass: 76kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

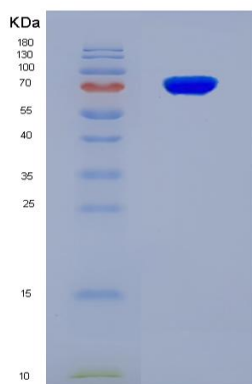


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.