

Recombinant Ras-related protein Rab-27A (Rab27a)

Catalog No.: TP11342

100µg

Sequence Information

Species: Mouse

Gene ID:11891

Swiss Prot:Q9ERI2

Synonyms:RB27A

Residues:Met1-Cys221

MSDGDYDYLIKFLALGDSGVGKTSVLYQYTDGKFNSKFITTVGIDFREKRVVYR
ANGPDGAVGRGQRIHLQLWDTAGQERFRSLTTAFFRDAMGFLLLFDLTNEQSFL
NVRNWISQLQMHAYCENPDIVLCGNKSDLEDQRAVKEEEEARELAEKYGIPYFET
SAANGTNISHAIEMLLDLIMKRMERCVDKSWIPEGVVRNNGHTSADQLSEEKEK
GLCGC

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His and SUMO Tag

Subcellular Location: Secreted.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.0

Predicted Molecular Mass: 41.0kDa

Accurate Molecular Mass: 41kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

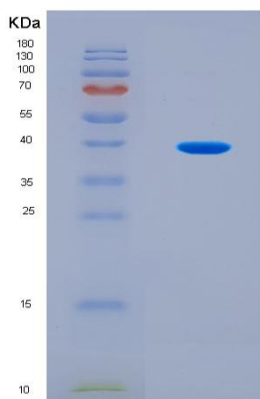


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.