

Recombinant Zinc finger and BTB domain-containing protein 20 (Zbtb20)

Catalog No.: TP10419

100µg

Sequence Information

Species: Mouse

Gene ID:56490

Swiss Prot:Q8K0L9

Synonyms:ZBT20

Residues:Met1-Gly741

MLERKKPKTAENQKASEENEITQPGGSSAKPALPCLNFEAVLSPAPALIHSTHS
LTNSHAHTGSSDCDISCKGMTERIHSINLHNFNSVLETLNEQRNRGHFCDVTV
RIHGSMRAHRCVLAAGSPFFQDKLLGYSDIEIPSVSVQSVQKLIDFMYSGV
LRVSQSEALQILTAASILQIKTVIDECTRIVSQNVGDVFPGIQDSGQDTPRGTP
ESGTSQGSSDTESGYLQSHQHSVDRIYSALYACSMQNGSGERSFYSGAVVSHH
ETALGLPRDHHMEDPSWITRIHERSQQMERYLSTTPEPTHCRKQPRPVRIQTLV
GNIHIKQEMEDDYDYYGQQRVQILERNESECTEDTDQAEGTESEPKGESFDSG
VSSSIGTEPDSEVQFGAAAPRDGQAEPAPQEAEEAPAESSAQPNQLEPGASS
PERSNESEMDNTVITVSNSSDKGVLQQPSVNTSIGQPLPSTQLYLRTETLTSN
LRMPLTLTNTQVIGTAGNTYLPALFTTQPAGSGPKPFLFSLPQPLTGQQTQFV
TVSQPGLSTFTAQLPAPQLASSAGHSTASGQGDKKPYECTLCNKFTAKQNYV
KHMVFHTGEKPHQCSICWRSFSLKDYLKHMVTHTGVRAYQCSICNKRFTQKSS
LNVHMLRHGEKSYECYICKKKFSHKTLLEHVALHSASNGTPPAGTPPGARAG
PPGVVACTEGTTYVCSVCPAKFDQIEQFNDHMRMHVSDG

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His and SUMO Tag

Subcellular Location: Nucleus.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.5

Predicted Molecular Mass: 97.1kDa

Accurate Molecular Mass: 97kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

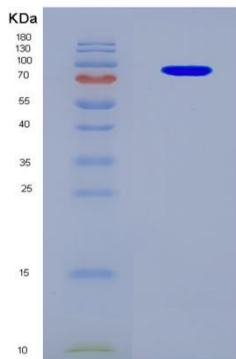


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.