

Recombinant Cytochrome P450 2A13 (CYP2A13)

Catalog No.: TP11398

100µg

Sequence Information

Species: Human

Gene ID:1553

Swiss Prot:Q16696

Synonyms:CP2AD

Residues:Gly31-Arg494

GKLPPGPTPLPFIGNYLQLNTEQMYNSLMKISERYGPVFTIHLGPRRWVVLGCH
DAVKEALVDQAEFFSGRGEQATFDWLFKGYGVAFSNGERAKQLRRFSIATLRGF
GVGKRGIEERIQQEAGFLIDALRGTHGANIDPTFFLSRTVSNVISSIVFGDRFD
YEDKEFLSLLRMMLGSFQFTATSTGQLYEMFSSVMKHLPGPQQQAFKELQGLED
FIAKKVEHNQRTLDPNSPRDFIDSFLIRMQEEKNPNTFYLKNLVMTTNLNFF
AGTETVSTTLRYGFLLLKMHPVEAKVHEEIDRVIGKNRQPKFEDRAKMPYTEA
VIHEIQRFGDMLPMGLAHRVNKDKFRDFFLPKGTEVFPMLGSVLRDPRFFSNP
RDFNPQHFLDKKGQFKKSDAFVPFSIGKRYCFGEGLARMELFLFFTITMQNFRF
KSPQSPKDIDVSPKHVGFATIPRNYTMSFLPR

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His and SUMO Tag

Subcellular Location: Endoplasmic reticulum membrane.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 9.5

Predicted Molecular Mass: 69.3kDa

Accurate Molecular Mass: 69kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

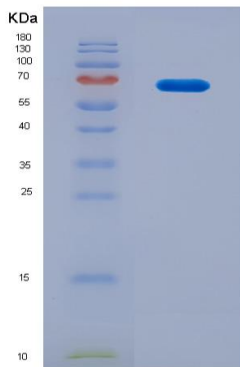


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.