



## Recombinant Human ALB Protein

Catalog No.: TP01119

50µg

### Sequence Information

**Species:** Human

**Gene ID:**85302

**Swiss Prot:**P02768

**Synonyms:**Serum albumin , PRO0883, PRO0903, PRO1341

**Residues:**Arg19-Leu609

RGVFRDDAHKSEVAHRFKDLGEENFKALVLI AFAQYLQQCPFEDHVKL VNEVTEFAKTCVADESAENC DKSLHTLFGDKLCTVATLRETYGEMADCCAKQEP  
ERNECFLQHKDDNP NLPRLVRPEVDVMCTAFHDNEETFLKKYLYEIARRHPYFYAPELLFFAKRYKAAFTECCQAADKAACLLPKLDEL RDEGKASSAKQRL  
KCASLQKFGERAFAKAWAVARLSQRFPKAEFAEVSKLVTDLT KVHTECCHGDLLECADDRADLAKYICENQDISSKLKECCEKPLLEKSHCIAEVENDEMPA  
DLPSLAADFVESKDVCKNYAEAKDVFLGMFLYEYARRHPDYSVVL LRLAKTYETTLEKCCAAADPHECYAKVFDEFKPLVEEPQNLIKQNC ELFQGEYK  
FQNALLVRYTKKVPQVSTPTLVEVSRNLGKVGSKCKHPEAKRMPCAEDYLSVVLNQLCVLHEKTPVSDRVTKCTESLVNRRPCFSALEVDETYVPKEFNA  
ETFTFHADICTLSEKERQIKKQTALVELVKHKPKATKEQLKAVMDDFAAFVEKCKADDKETCFAEEGKKLVAASQAALGL

### Product Information

**Source:** Eukaryotic expression.

**Host:** *Insect cell*

**Tags:** NoneTag

**Subcellular Location:** Extracellular matrix.

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 5% Trehalose .

**Original Concentration:** 1000µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 5.9

**Predicted Molecular Mass:**69kDa

**Accurate Molecular Mass:** 69kDa as determined by SDS-PAGE reducing conditions.

### [ USAGE ]

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

### [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

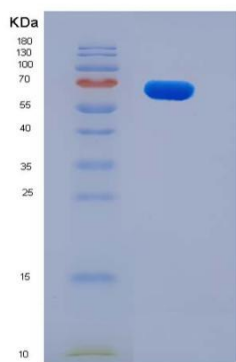


Figure 1. SDS-PAGE

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.