

## Recombinant Mitochondrial antiviral-signaling protein (MAVS)

Catalog No.: TP11458

100µg

### Sequence Information

**Species:** Human

**Gene ID:**57506

**Swiss Prot:**Q7Z434

**Synonyms:**MAVS

**Residues:**Pro2-Pro513

PFAEDKTYKYICRNFSNFCNVDVVEILPYLPCLTARDQDRLRATCTLSG NRDTL  
WHLFNTLQRRPGWVEYFIAALRGCELVDLADEVASVYQSYQPRTSDRPPDPLEP  
PSLPAERP GPPTPAAAH SIPYNSCREKEPSYMPVQETQAPESPGENSEQALQT  
LSPRAIPRNP DGGPLESSSDLAALSPLTSSGHQE QDTELGSTHTAGATSSLTPS  
RGPVSPSVSFQPLARSTPRASRLPGPTGSVVSTGT SFSSSSPGLASAGAAEGKQ  
GAESDQAEPIICSSGAEAPANSLPSKVPTT LMPVNTVALKVPANPASVSTVPSK  
LPTSSKPPGAVPSNAL TNPAPSKLPINSTRAGMVPSKVPTSMVLTKVSASTVPT  
DGSSRNEETPAAPT PAGATGGSSAWLDSSSEN RGLGSEL SKPGVLASQVDSPFS  
GCFEDLAISASTSLGMGPCHGPEENEYKSEGTFGIHVAENPSIQ LLEGNP GPPA  
DPDGGPRPQADRK FQEREVPCHRPSP

### Product Information

**Source:** Recombinant expression.

**Host:** *E.coli*

**Tags:** N-terminal His and SUMO Tag

**Subcellular Location:** Mitochondrion outer membrane.

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 4.8

**Predicted Molecular Mass:** 69.4kDa

**Accurate Molecular Mass:** 69kDa as determined by SDS-PAGE reducing conditions.

### **[ USAGE ]**

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

## **[ STORAGE AND STABILITY ]**


**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## **[ IDENTIFICATION ]**

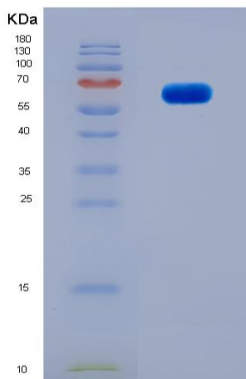

**Figure 1. SDS-PAGE**

## **[ IMPORTANT NOTE ]**

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.