

Recombinant Mouse Interleukin 21 Receptor (IL21R)

Catalog No.: TP11537

50µg

Sequence Information

Species: Mouse

Gene ID:60504

Swiss Prot:Q9JHX3

Synonyms: CD360 antigen; CD360; IL-21 R;

IL-21 receptor; IL21R; IL-21R;

interleukin 21 receptor; interleukin-21

receptor; MGC10967; NILR; Novel

interleukin receptor

Residues:Cys20 - Pro236

CLDLTCYTDYLWTITCVLETRSPNPSILSLTWQDEYEELQDQETFCSLHRSGHN

TTHIWYTCMRLSQFLSDEVFIVNVTQSGNNSQECGSFVLAESIKPAPPLNVT

VAFSGRYDISWDSAYDEPSNYVLRGKLQYELQYRNLRDPYAVRPVTKLISVDSR

NVSLLP EEFHKDSSYQLQVRAAPQPGTSFRGTWSEWSDPVIFQTQAGEPEAGWD

P

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: C-terminal His and FC Tag

Subcellular Location: Membrane.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 1000µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 4.4

Predicted Molecular Mass: 51.5kDa

Accurate Molecular Mass: 70kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

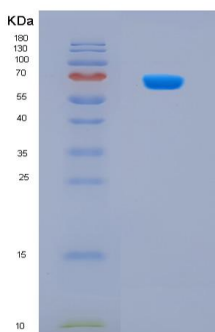


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.