

Recombinant Cytochrome P450 3A4 (CYP3A4)

Catalog No.: TP10606

100µg

Sequence Information

Species: Human

Gene ID:1576

Swiss Prot:P08684

Synonyms:CP3A4

Residues:Ala2-Ala503

ALIPDLAMETWLLLA VSLVLLYLYGTHSHGLFKKLGI PGTPLPFLGNILSYHK
GFCMFDMECHKKYGKVG FYDGGQPVLAITDPDMIKT VLVKECYSVFTNRRPFG
PVGFMKSAISIAEDEEWKRLRSLLSPTFTSGK LKEMVPIIAQYGDVLVRNLRRE
AETGKPVTLKDVFGAYSMDVITSTSGV NIDSLNNPQDPFVENTKKLLRFDFLD
PFFLSITVPFLIPILEV LNICVFPREVTNFLRKSVKRMKESRLED TQKHRVDF
LQLMIDSQNSKETESHKALSDLELVAQSIIFIFAGYETTSSVLSFIMYELATHP
DVQQKLQEEIDAVLPNKAPPTYD TVLQMEYLD MVVNETLRLFP IAMRLERVCKK
DVEINGMFIPKGVVVMIPSYALHRDPKYWTEPEKFLPERFSKKNKDNIDPYIYT
PFGSGPRNCIGMRFALNMNKLALIRVLQNFSFKPKETQIPLKLSLGGLLQPEK
PVVLKVESRDGTVSGA

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His and SUMO and AVI Tag

Subcellular Location: Endoplasmic reticulum membrane.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.2

Predicted Molecular Mass: 74.2kDa

Accurate Molecular Mass: 74kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

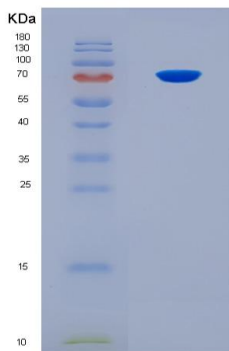


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.