

Recombinant Cleavage and polyadenylation specificity factor subunit 2 (CPSF2)

Catalog No.: TP11558

100µg

Sequence Information

Species: Human

Gene ID:53981

Swiss Prot:Q9P2I0

Synonyms:HDAC1; GON10; HD1; RPD3;
RPD3L1

Residues:Met1-Val782

MTSIIKLTTLSGVQEEALCYLLQVDEFRLFLLDCGWDEHFSMDIIDSLRKHHVHQ
IDAVLLSHPDPLHLGALPYAVGKLGLENCAIYATIPVYKMGQMFMYDLYQSRHNT
EDFTLFTLDDVDAAFDKIQQLKFSQIVNLKGKGHGLSITPLPAGHMIGGTIWKI
VKDGEIEIVYAVDFNHKREIHLNGCSLEMLSRPSLLITDSFNATYVQPRRKQRD
EQLLTNVLETLRGDGNVLIADVTAAGRVLLEAQLLDQIWRTKDAGLGVYSLALLN
NVSYNVVEFSKQVWMSDKLMRCFEDKRNNPFQFRHLSLCHGLSDLARVPSPK
VVLASQPDLECGFSRDLFIQWCQDPKNSIILTYRTTPGTLARFLIDNPSEKITE
IELRKRVKLEGKELEEYLEKEKLKKEAAKKLEQSKEADIDSSDESDEEDIDQP
SAHKTKHDLMMKGEGSRKGSFFKQAKKSYPMFPAPEERIKWDEYGEIIKPEDFL
VPELQATEEEKSKLESGLTNGDEPMDQDLSDVPTKCISTTESIEIKARVTYIDY
EGRSDGDSIKKIINQMKPRLIIVHGPPEASQDLAECCRAFGGKDIKVMPKLH
ETVDATSETHIYQVRLKDSLVSLLQFCKAKDAELAWIDGVLDMRVSKVDTGVIL
EEGELKDDGEDSEMQVEAPSDDSVIAQQKAMKSLFGDDEKETGESEIIPITLPE
LPPHEVPGHQSVFMNEPRLSDFKQVLLREGIQAEFVGGVLVCNNQVAVRRRTETG
RIGLEGCLCQDFYRIRDLLYEQYAIIV

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His and SUMO Tag

Subcellular Location: Nucleus.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 4.7

Predicted Molecular Mass: 104.5kDa

Accurate Molecular Mass: 105kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

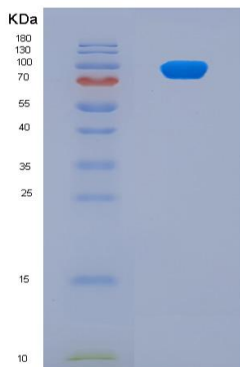


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.