

Recombinant Human Albumin (ALB)

Catalog No.: TP10680

50µg

Sequence Information

Species: Human

Gene ID:213

Swiss Prot:P02768

Synonyms: ALBU; Serum albumin

Residues:Asp25-Leu609

DAHKSEVAHRFKDLGEENFKALVLIAFAQYLQQCPFEDHVKLNVTEFAKTCV
ADESAENCDSLHTLFQDKLCTVATLRETYGEMADCCAKQEPERNECF LQHKDD
NPNLPRLVRPEVDMCTAFHDNEETFLKKYLYE IARRHPYFYAPELLFFAKRYK
AAFTGCCQAADKAACLLPKLDEL RDEGKASSAKQRLKCASLQKFGERAFAKAWAV
ARLSQRFPKAEFAEVSKLVTDLT KVHTECCHGDLLECADDRADLAKYICENQDS
ISSKLKECCEKPLLEKSHCIAEVENDEMPADLPSLAADFVESKDVCKNYAEAKD
VFLGMFLY EYARRHPDYSVVL LRLAKTYETTLEKCCAAADPHECYAKVFDEFK
PLVEEPQNLIKQNCLEQLGEYKFQNALLVRYTKKVPQVSTPTLVEVSRNLGK
VGSCKCKHPEAKRMPCAEDYLSVVLNQ LCVLHEKTPVSDRVT KCCTESLVNRRP
CFSALEVDETYVPKEFNAETFTFHADICTLSEKERQIKKQTALVELVKHKPKAT
KEQLKAVMDDFAAFVEKCKADDKETCFAEEGKKLVAA SQAALGL

Product Information

Source: Eukaryotic expression.

Host: *Mammalian cells*

Tags: C-terminal Fc-Tag

Subcellular Location: Secreted.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 1000µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.8

Predicted Molecular Mass: 76.5kDa

Accurate Molecular Mass: 77kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

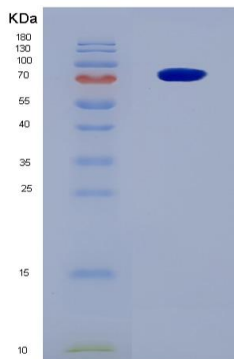


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.