



## Recombinant Human TMPRSS2 Protein

Catalog No.: TP11665 100µg

### Sequence Information

**Species:** Human

**Gene ID:** 7113

**Swiss Prot:** O15393

**Synonyms:** Serine protease 10, PRSS10, TMPRSS2, Transmembrane protease serine 2

**Residues:** Trp106-Gly492

WKFMGSKCSNSGIECDSSGTCINPSNWCDGVSHCPGGEDENRCVRLYGPNFILQVYSSQRKSWHPVCQDDWENYGRAACRDMGYKNNFYSSQGIVDD  
SGSTSFMKLNTSAGNVDIYKKLYHSDACSSKAVVSLRCIACGVNLNSSRQSRIVGGESALPGAWPWQVSLHVQNVHVCGGSIITPEWIVTAAHCVEKP  
LNNPWHWTAFAGILRQSFMYGAGYQVEKVISHPNYDSKTKNNDIALMKLQKPLTFNDLVKPVCLPNPGMMLQPEQLCWISGWGATEEKGKTSEVLNA  
AKVLLIETQRCNSRYVDNLITPAMICAGFLQGNVDSQGDGGPLVTSKNNIWWLIGDTSWGS GCAKAYRPGVYGNVMVFTDWIYRQMRADG

### Product Information

**Source:** Eukaryotic expression.

**Host:** Mammalian cells

**Tags:** C-terminal Fc-Tag

**Subcellular Location:** / .

**Purity:** >95%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, 5% Trehalose.

**Original Concentration:** 1000µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 8.12

**Predicted Molecular Mass:** 70.0kDa

**Accurate Molecular Mass:** 71kDa as determined by SDS-PAGE reducing conditions.

### **[ USAGE ]**

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

### **[ STORAGE AND STABILITY ]**

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

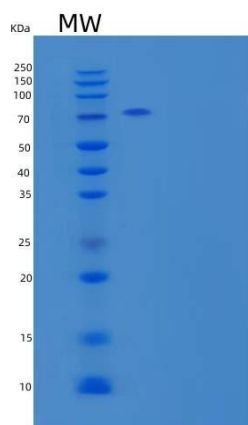


Figure 1. SDS-PAGE

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.