

Recombinant DNA ligase (ligA)

Catalog No.: TP10778

100µg

Sequence Information

Species: ECOLI

Gene ID:946885

Swiss Prot:P15042

Synonyms:DNLJ

Residues:Met1-Ser671

MESIEQQLTELRTTLRHHEYLYHVMDAPEIPDAEYDRLMRELRELETKHPELIT
PDSPTQRVGAAPLAAFSQIRHEVPMLSLDNVFEESFLAFNKRVDRLKNNKEV
TWCCCLKLDGLAVSILYENGLVLSAATRGDGTGEDITSNVRTIRAIPCLKHGE
NIPARLEVRGEVFLPQAGFEKINEDARRTGGKVFANPRNAAAGSLRQLDPRITA
KRPLTFFCYGVGVLEGGELPDTHLGRLLQFKKWGLPVSDRVTLCESAEEVLAFY
HKVEEDRPTLGFIDIDGVVIVNSLAQQEQLGFVARAPRWAVAFKFFPAQEOMTFV
RDVEFQVGRGTGAITPVARLEPVHVAGVLVSNATLHNADEIERLGLRIGDKVVIR
RAGDVIPQVNVVLSERPEDTREVFPTHCPVCGSDVERVEGEAVARCTGG LIC
GAQRKESLKHVFSRRAMDVDGMGDKIIDQLVEKEYVHTPADLFKLTAGKLTGLE
RMGPKSAQNVVNALEKAKETTFARFLYALGIREVGEATAAGLAAYFGTLEALEA
ASIEELQKVPDVGIVVASHVHNFFAEESNRNVESELLAEGVHWPAPIVINAE EI
DSPFAGKTVVLTGSLSQMSRDDAKARLVELGAKVAGSVSKKTDLVIAGEAAGSK
LAKAQELGIEVIDEAEMLRLLGS

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His Tag

Subcellular Location: /.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.2

Predicted Molecular Mass: 77.0kDa

Accurate Molecular Mass: 77kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

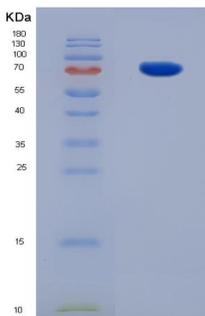


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.